

PAPETTE® NON-STERILE

For Cervical Cytology utilizing a Liquid-based cytology (LBC) method



Read all Safety information and instructions listed below before using this product!

Caution: U.S federal law restricts this device to sale by or the order of a physician.


QTY 100 Wallach Papettes® per bag (500/box)

REF LGM Catalog # 90-500

Application: The *PAPETTE*® is the most gentle, yet effective means of simultaneously taking both ectocervical and endocervical cell samples.

1. Contact the cervix with the *PAPETTE*® and insert the central bristles into the cervical canal deep enough to allow the shorter bristles to fully contact the ectocervix. Be sure not to force the shorter bristles into the cervical os to prevent unintended separation of the brush from the handle.
2. While maintaining gentle pressure in the direction of the cervix, rotate the brush two or more times in either direction.
3. Open the Preservative vial, detach the brush head by gently pushing it into the vial. Place the cap back on the vial and hand tighten. Mark the vial label with pertinent patient information. Forward the vial to the laboratory.
4. Dispose of *PAPETTE*® in accordance with local Medical Hazardous Waste practices.

Be sure not to force the shorter bristles into the cervical os to prevent unintended separation of the brush from the handle.

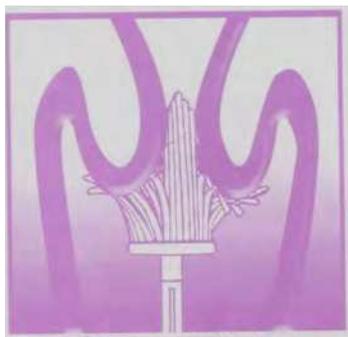

PAPETTE® U.S Patent #4,700,713 Latex Free

Distributed By:
LGM International Inc.
3030 Venture Lane, Suite 106
Melbourne ,Florida 32934
Telephone: (321) 254-0480
Fax: (321) 254-0481

Papette® is a registered trademark of :
Wallach Surgical Devices
Orange CT, USA.

Authorized European Representative:
Leisegang Feinmechanik GmbH
Leibnizstrasse 32
D-10625, Berlin GERMANY



PAPETTE® NON-STERILE

For Cervical Cytology utilizing a Liquid-based cytology (LBC) method



Read all Safety information and instructions listed below before using this product!

Caution: U.S federal law restricts this device to sale by or the order of a physician.


QTY 100 Wallach Papettes® per bag (500/box)

REF LGM Catalog # 90-500

Application: The *PAPETTE*® is the most gentle, yet effective means of simultaneously taking both ectocervical and endocervical cell samples.

1. Contact the cervix with the *PAPETTE*® and insert the central bristles into the cervical canal deep enough to allow the shorter bristles to fully contact the ectocervix. Be sure not to force the shorter bristles into the cervical os to prevent unintended separation of the brush from the handle.
2. While maintaining gentle pressure in the direction of the cervix, rotate the brush two or more times in either direction.
3. Open the Preservative vial, detach the brush head by gently pushing it into the vial. Place the cap back on the vial and hand tighten. Mark the vial label with pertinent patient information. Forward the vial to the laboratory.
4. Dispose of *PAPETTE*® in accordance with local Medical Hazardous Waste practices.

Be sure not to force the shorter bristles into the cervical os to prevent unintended separation of the brush from the handle.

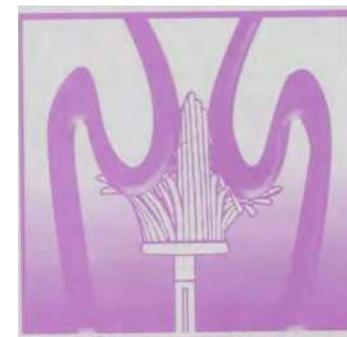

PAPETTE® U.S Patent #4,700,713 Latex Free

Distributed By:
LGM International Inc.
3030 Venture Lane, Suite 106
Melbourne ,Florida 32934
Telephone: (321) 254-0480
Fax: (321) 254-0481

Papette® is a registered trademark of :
Wallach Surgical Devices
Orange CT, USA.

Authorized European Representative:
Leisegang Feinmechanik GmbH
Leibnizstrasse 32
D-10625, Berlin GERMANY

